PROGENT

Intensive cleaner

INDICATIONS

Protein remover and disinfecting solution for all rigid gas permeable lenses.

DIRECTIONS

- 1 Clean your lenses thoroughly with regular cleaning solution.
Place the lenses into vial caps in its appropriate right or left holder : R = right, L = left, for intensive cleaning.
- 2 Open both doses carefully, by holding upright. Gently and slowly twist the tip, without touching the opening or squeezing the dose.
- 3 Insert the tip of dose inside the vial and fill slowly, first with A solution, then B solution, taking care not to spill any solution.
- 4 Close the vial cap securely gently shake it a few times.
- 5 Soak the lenses in Progrent for 30 minutes maximum.
- 6 Remove the vial cap, empty its content. Rinse the lenses in the holder thoroughly with the regular rinsing solution.
- 7 Then clean the lenses thoroughly with the regular solution and rinse them. The lenses may then be worn or stored in the storage case.
- 8 Rinse the vial and its cap thoroughly with the rinsing and regular and let them air dry.

WARNINGS

- Do not apply Progrent to the eyes, the nose and the mouth. In case of misuse, rinse immediately your eyes thoroughly with water and consult your eye care practitioner.
- It is essential to rinse the lenses abundantly with a rinsing and regular solution after soaking them in Progrent solution.
- KEEP OUT OF REACH OF CHILDREN

PRECAUTIONS

- Do not soak lenses in Progrent for over 30 minutes. Lenses might discolour. However this will not affect their quality.
- Important: do not swallow.
- To be used weekly for daily use (or on removing extended wear lenses).
- Do not use Progrent A and B solutions separately.
- Single use doses. Discard the solution and the doses after each use.
- Store at room temperature (15-25°C).
- Avoid contact with the skin and do not spill Progrent on your clothes. It may bleach them.
- Do not use Progrent A and B after expiry date printed on the single dose.
- Store out of direct sunlight.
- Do not use Progrent A and B doses if they are damaged.
- Consult your eye care professional before using any medication with contact lenses or before changing solution.

CONTRA-INDICATIONS

- Do not use for soft lenses.
- Do not use if allergic to one of the ingredients.

ADVERSE REACTIONS

If you experience persistent discomfort or irritation, remove your lenses, stop using the product and consult your eye care professional.
Any serious incident that has occurred in relation to the device should be reported to the manufacturer and the competent authority of the Member State in which the user and/or patient is established.

CONTACT LENS DISINFECTION


Progrent acts as a disinfectant for microorganisms (bacteria, moulds, viruses, amoebae) in just 5 minutes of soaking.

PACKAGING

Progrent is available in boxes of 5 pairs of Progrent A and B with one Progrent vial and starter kit.

COMPOSITION

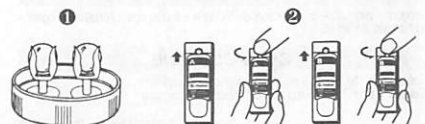Solution A : NaClO 20.3 mg, purified water ad 5ml/dose.
Solution B : KBr 30.75 mg, purified water ad 5ml/dose.

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